2019-nCoV Antigentest (Lateral-Flow-Methode)

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**Wichtig**

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# PROCEDURE CARD

**Important**

- Nasal swab specimens from anyone under the age of 18 should be carried out under adult supervision. People over the age of 75 should be assisted during the swabbing procedure.
- The test takes about 5 minutes to complete, and the test results must be read after 15 minutes; to ensure test accuracy, do not read results after 20 minutes.
- It is important that you read the test instructions for use carefully and follow the test procedures to get an accurate result.
- It is recommended to perform the test at room temperature (15-30°C).
- Please wash or disinfect your hands carefully before and after performing the test.
- Please take the necessary precautions when testing other people (e.g. face mask, gloves).

**Product Components**

1. 2019-nCoV Antigen Test Cassette
2. Pre-filled Extraction Buffer
3. Disposable Sterile Swab
4. Biohazard Specimen Bag
5. Instruction for Use
6. Procedure Card


Other require items (not in the self-test included)

- Clock or timer


Illustration exemplary


**Test Procedure**

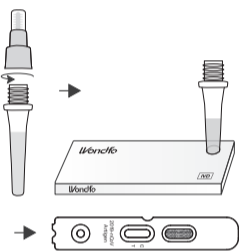**Step 1**  
Take out the Pre-filled Extraction Buffer Tube, unscrew the lid and place the tube in the tube rack (The tube rack is on or inside the test kit).  
Take out the Test Cassette from the foil pouch and lay it out flat.

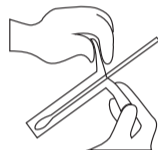
**Step 2**  
Remove the swab from the container, being careful not to touch the soft end, which is the absorbent tip.

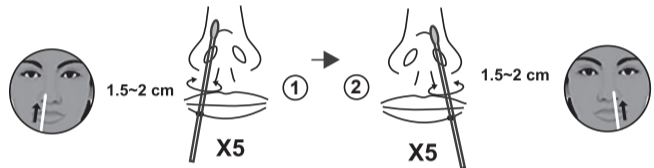
**Step 3**  
Carefully insert the entire absorbent tip of the swab into a nostril (about 1.5~2 cm). Take samples from the nasal wall by firmly rotating the swab in a circular path five times against the nasal wall. Slowly remove the swab from the nostril. (This step should take approximately 15 seconds to ensure the collection of mucous and cells). Repeat the sampling as described above for the other nostril with the same swab.


**Note:** Simply twirling the swab against one part of the inside of the nose or leaving the swab in the nose for 15 seconds is not a proper technique and may result in an insufficient sample.

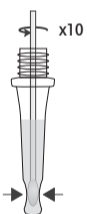
**Step 4**  
Insert the swab into the Pre-filled Extraction Buffer Tube and immerse the entire tip of swab into the Extraction Buffer. Rotate about 10 times.

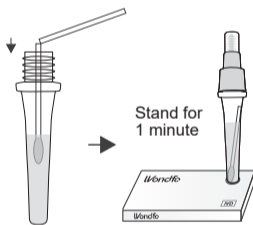**Step 5**  
Gently take out the swab, breaking off the tip of swab from the first break-point. Leave it in the tube, rescrew the cap and place the tube into the tube rack for 1 minute.

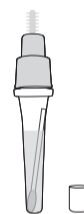**Step 6**  
Unscrew the small cap at the top of the Pre-filled Extraction Buffer Tube.

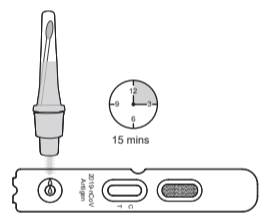**Step 7**  
Lay the Test Cassette flat and add 4 drops of processed specimen into the sample well. Wait for 15 minutes to get the result. **Do not read results after 20 minutes.**


**Step 8**  
After the test is completed, place the Test Cassette in the Biohazard Specimen Bag and dispose all test kit materials via biohazard waste disposal protocol.


**Step 9**  
Re-apply hand sanitizer


**Step 10**  
Result interpretation

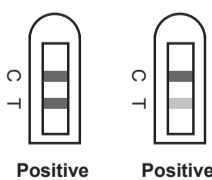Colored bands appear at both test line (T) and control line (C). This indicates a positive result for the 2019-nCoV antigen in the specimen. **It does not matter if one of the lines that make up the test line (T) is lighter or darker than other; the result is "Positive".**

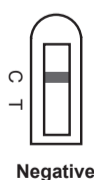Colored band appears at control line (C) only. This indicates that the concentration of the 2019-nCoV antigen is zero or below the detection threshold of the test.

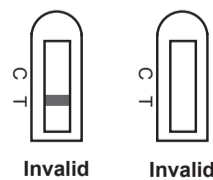No visible colored band appears at control line after performing the test. The instructions may have not been followed correctly or the test may have deteriorated. It is recommended to repeat the sampling and testing procedures.